

OxyHome® Quick Reference Guide

1. Using the Pause Button

Push Pause when taking a short break from the concentrator. This feature stops oxygen delivery and continues to cool the compressor to improve longevity. Push pause again to resume operation.

2. Alarm On/Off Button

Push and hold for 10 seconds to mute or unmute the alarm.

Please note: If the alarm is muted, it will not be audible until unmuted or until the machine is turned off and back on.

3. Timer Button

Push and hold for 10 seconds to activate timer function. The timer can be set up to 300 minutes in 10 minute increments.


4. Checking the Purity

Press and hold the "Pause" button for 5 seconds to display Purity.

5. Indicator Lights

There are 4 indicator lights on the device:

1. Low outlet pressure indicator (Low Pressure): Outlet pressure is less than 2.2 psi, the indicator light will flash yellow every two seconds and an alarm will sound.
2. Low oxygen purity indicator (Low Conc.): If the oxygen concentration drops below 82% (V/V) the indicator light will flash yellow, and the alarm will sound.
3. Power loss indicator (Power Outage): When there is a power outage, the indicator light will flash yellow every 2 seconds and the alarm will sound.
4. Power status indicator: When the concentrator is powered on, the indicator light will remain solid green if the concentrator is working normally.

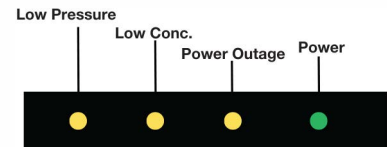

6. Alarms

Item No.	Alarm type	Priority	Alarm Signals / Indicator light			Auditory Yes/No	Alarm condition delay/s	Alarm signal generation delay/s
			Color	Flashing frequency	Vacancy ratio			
1	Power supply failure alarm	Med	Yellow	Flashing every two seconds	50%	Yes	5s	5s
2	Start-up fail alarm	Med	Yellow	Flashing every two seconds	50%	Yes	5s	5s
3	Low oxygen concentration ALARM (non-startup)	Med	Yellow	Flashing every two seconds	50%	Yes	5 s	5s
4	Pressure failure alarm	Med	Yellow	Flashing every two seconds	50%	Yes	5 s	5s
5	No flow alarm	Med	On screen	Flashing every two seconds	N/A	Yes	5 s	5s
6	Low flow alarm	Med	On screen	Flashing every two seconds	N/A	Yes	10 s	5s

7. Replacement Filters

Part Number	Description
1170-5401	Replacement foam cabinet filter for OxyHome 5 Liter Stationary Concentrator
1170-5402	Replacement HEPA/Compressor inlet filter for OxyHome 5 Liter Stationary Conc.

