

**FREESTYLE® COMFORT®**

**A NEW STANDARD IN  
PORTABLE OXYGEN.**



**Smart O<sub>2</sub> delivery.  
Unmatched comfort.  
Durable design.**

**CPAPX**

  
**CAIRE®**  
*O<sub>2</sub> is what we do.*

# INNOVATION MEETS OXYGEN THERAPY.

Expand your patients' comfort zone with this innovative and powerful device. Its advanced clinical performance makes it easy for patients to stay saturated, active and healthy.

## Clinically-Effective Oxygen Delivery



### Smart technology to ensure saturation

- autoSAT®** Automatic titration to patient's respiratory rate
- autoDOSE®** Automatic pulse delivery when no breath detected
- UltraSense®** Sensitive breath detection
- myCAIRE®** Telehealth solution



### Flexible flow rates to adjust to changing needs

With 5 settings and up to 1050 mL of oxygen per minute, the FreeStyle Comfort offers flexibility in treatment options and covers broad range of prescriptions.

## Unmatched Comfort and Convenience



### Compact and Lightweight

At only 25.4 cm tall and 2.3 kg, it is similar in size to a book and weighs less than a purse.



### Quiet

At 39.9 dBA the sound level is slightly more than a whisper and less than a refrigerator.\*\*



### Ergonomically Designed

Body-skimming, curved silhouette reduces swinging and bouncing.



### Easy-to-Read Display

Large fonts, multi-color LCD display and glow-in-the-dark keypad.



### Lasts All Day

The 8-cell battery pack lasts up to 8 hours and the 16-cell battery pack lasts up to 16 hours.\*



### Enables Airline Travel

Meets FAA requirements for commercial airline travel.

## Proven Quality and Support



### Industry-Leading Reliability

Greater than 99% reliability rate out of the box.



### Reduced Maintenance

Robust design features superior sieve bed technology reducing required maintenance.



### Comprehensive Long-Term Warranty

CAIRE is the first manufacturer to offer a standard 3-year warranty on the device AND sieve beds.



# PATIENT SATISFACTION BUILT-IN.


Smart features  
to keep patients  
saturated

Ergonomic curved design  
to let patients move freely

Large, easy-to-read,  
color display

Small, lightweight (2.3 kg.)

Long lasting battery life –  
enough to let patients stay out all day

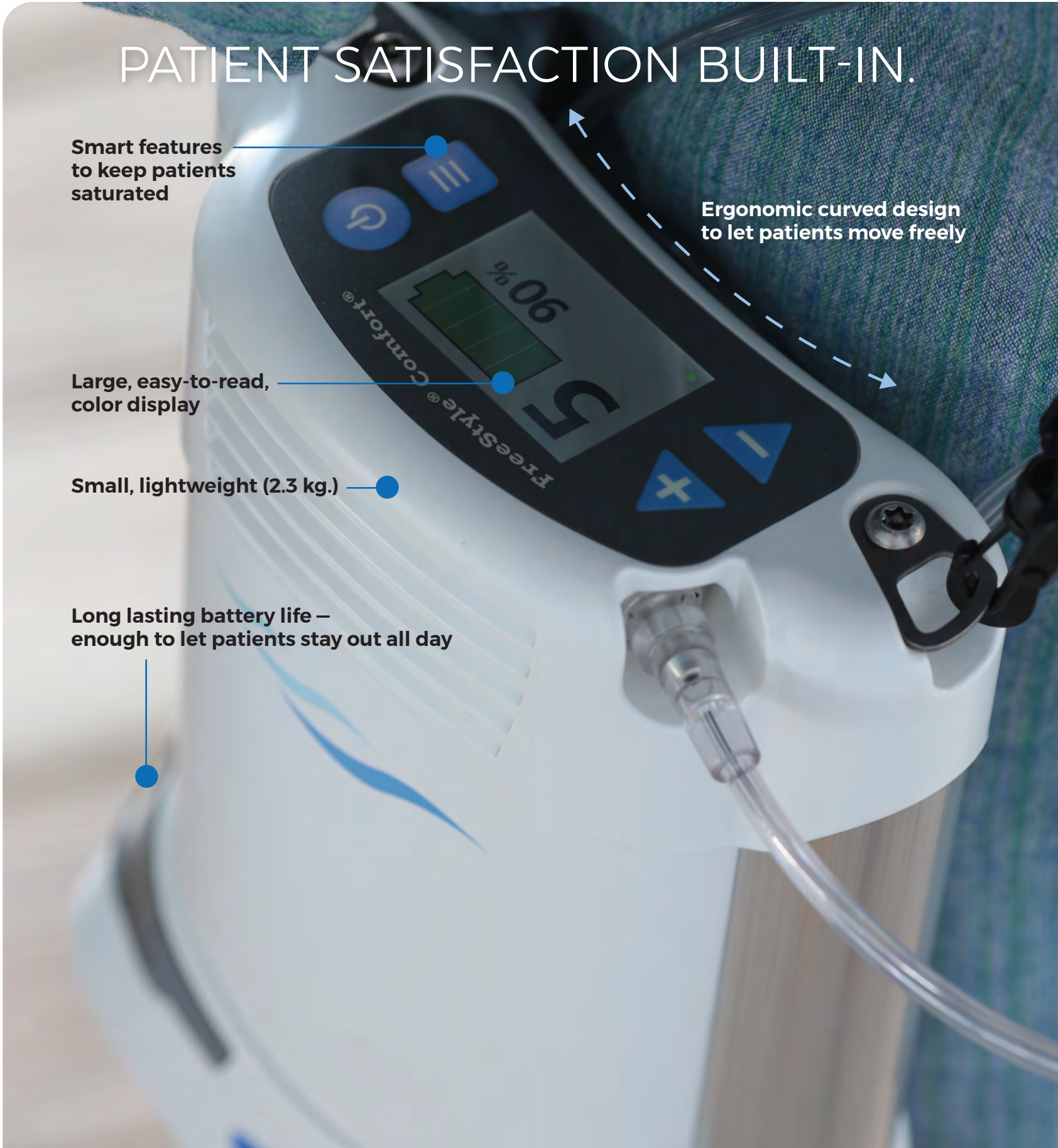

# FreeStyle® Comfort®

The perfect fit



## Key features and benefits



### Optimal Comfort

Let your patients experience the comfort of the body-skimming, curved silhouette design. The FreeStyle Comfort portable oxygen concentrator is lightweight and compact.



### Smart O2 delivery

Unique clinical features including autoSAT®, UltraSense®, autoDOSE®, and myCAIRE™ help to ensure your patients will get the oxygen they need, when they need it.



### Proper O2 Saturation

With 5 settings and up to 1050 mL of oxygen per minute, the FreeStyle Comfort serves a wide range of patient needs. The unit delivers up to 95% pure oxygen to further ensure proper saturation.



### Proven Quality

Robust design features, including the sieve bed technology and critical components, allow for minimum maintenance and long life. CAIRE offers one of the best warranties available, with a 3-year warranty on the device and sieve beds.

## Specifications

### FLOW RATE

Pulse Dose 1-5  
O2 Output Settings:  
Setting 1-210 mL/min  
Setting 2-420 mL/min  
Setting 3-630 mL/min  
Setting 4-840 mL/min  
Setting 5-1050 mL/min

### DIMENSIONS

10.0 x 7.3 x 3.1 in  
(25.4 x 18.5 x 7.9 cm)  
8-cell battery pack  
  
11.0 x 7.3 x 3.1 in  
(27.9 x 18.5 x 7.9 cm)  
16-cell battery pack

### WEIGHT

5 lbs (2.3 kg)  
8-cell battery pack  
  
6 lbs (2.7 kg)  
16-cell battery pack

### SOUND LEVEL

39.9 dB(A) @ Setting 2

### TRIGGER SENSITIVITY

>-0.10 cm/H2O  
(adjustable)

### BATTERY DURATION

8-cell battery  
up to 8 hours  
  
16-cell battery  
up to 16 hours

### POWER

AC Power: 100-240  
VAC, 50-60Hz  
DC Power: 11-18 VDC  
(10 max amp)  
Battery: Lithium Ion

### WARRANTY

3 years  
(Includes sieve beds)  
  
1 year on battery and  
power supplies

