

PARI Trek® S Portable Aerosol System

	Part #	NDC #	HCPCS #
Standard	047F45-LCS	83490-0470-06	E0570, A7005
Deluxe	047F35-LCS	83490-0470-05	E0571, A7005

All configurations include Wing Tip™ Tubing



Trek® S Premium Packaging



Deluxe Carrying Case



PARI LC® Sprint Reusable Nebulizer and backup reusable nebulizer



Timestrip® - 6 Month Reminder



Trek® S Compressor



AC Adapter & 12V DC Adapter



Trek® S Instructions & Instructional DVD




Li-Ion Rechargeable Battery
Deluxe Only


The Difference is in the Delivery

Trek® S vs. Other Portable Nebulizer Systems²

DRUG DELIVERED TO THE LUNGS

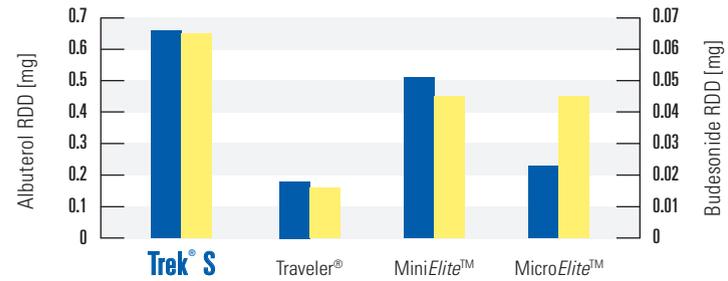

DRUG DELIVERED TO THE LUNGS PER MINUTE

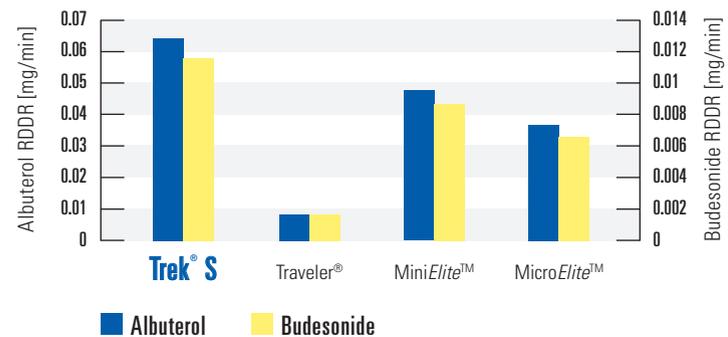

Performance with PARI LC® Sprint³

Total Output Rate (TOR) = **440 mg/min**
 Mass Median Diameter (MMD) = **3.8 µm**
 Percentage Respirable Fraction (%RF) = **65 %**
 Average Treatment Time (ATT) = **5-6 min**

Technical Specifications:

Maximum Pressure of Compressor = **35 PSI**
 Maximum Liter Flow of Compressor = **8.5 LPM**
 Nominal Operating Pressure⁴ = **14.5 PSI (1.0 bar)**
 Nominal Air Flow⁴ = **4.0 LPM**
 Weight = **0.8 lbs (1.2 lbs with Battery)**
 Compressor Warranty = **3 years**
 Battery Warranty = **6 months**
 Battery Life = **50 min**
 Time to Recharge = **2 hours**

Accessories:

Portable Power Pack = **047F25**
 PARI LC® Sprint = **023F35**
 Replacement Air Filters = **041E4851P12**

CPAPX



TREK® S

PORTABLE AEROSOL SYSTEM

From the leader in clinically proven aerosol delivery



PARI nebulizers are the most widely used nebulizers by pharmaceutical companies in clinical trials for new medications compared to any other brand.

WHY TREK® S?



“Fast treatments combined with ultimate portability allow me to enjoy an active lifestyle.”



“Because patients have the freedom to take treatments anywhere, I am confident they will adhere to their treatment plan.”



“I depend on the Trek® S for safe aerosol therapy no matter where we are - at home, school, or even the big game.”

LC® Sprint



Experience the freedom of portable aerosol therapy! The PARI Trek® S Portable Aerosol System is designed from PARI's long legacy of clinically proven products and incorporates fast, efficient aerosol treatments with portable technology. Perfect for those with active lifestyles, Trek® S allows patients the freedom to take treatments anywhere.

¹Lester, et al. Nebulizer Use and Maintenance by Cystic Fibrosis Patients: A Survey Study. *Respiratory Care* 2004; 49(12): 1504-1508.

²Measured with Sierra Cascade Impactor; Albuterol 2.5mg/3mL; Budesonide .25mg/2mL; Inspiratory Flow 30liters/min; Breathing pattern (500mLs, 15cycles/min, duty cycle 0.5)

³Measured with Malvern MasterSizer X at 50% relative humidity, 0.9% NaCl solution, inspiratory flow 20 liters/minute, continuous nebulization, 23° C, fill volume 2.5 ml.

⁴When used with a PARI LC® reusable nebulizer.

⁵All registered marks and trademarks belong to their respective companies.

Consistent Aerosol Delivery

EXCLUSIVE TO PARI, **TIMESTRIP®** IS AN EASY TO USE, VISUAL INDICATOR THAT CLEARLY INSTRUCTS PATIENTS TO REPLACE THEIR PARI REUSABLE NEBULIZERS EVERY 6 MONTHS. NO MORE GUESSWORK OR HASSLE WITH CALENDARS.



- Easy to use and activate
- Great for Patient Reminder programs
- Benefits of replacing PARI nebulizers:
 1. Consistent medication delivery
 2. Fast treatment times
 3. Reduces the risk for bacterial contamination¹